



Molift Quattro

The temporary hoist solution with excellent stability

A freestanding traverse gantry

Molift Quattro is a freestanding hoist unit that offers a (four post) traverse system, which is the best solution for temporary hoisting needs or when the facility is not built for fixed installations. The gantry is easy to install without the need to mount additional fixtures on walls or ceilings. It is ideally suited for temporary hoist applications in nursing or private homes.



Features



Supreme stability & manoeuvrability

Molift Quattro is a gantry system which has supreme stability and provides access in all lifting situations.



Easy to install

Molift Quattro is designed for quick installation by one or (preferably) two persons.



Charging of hoist motor

Molift hoist motors can easily be charged through the hand control or through the in-rail charging option.

Solves most overhead lifting situations

The traverse has a lightweight aluminum frame and four steel posts that provides excellent stability and covers the surface of the hoisting requirement. Molift Quattro solves all overhead lifting situations together with the Molift Air or Molift Nomad hoist motor and a sling.

With a lifting capacity of max 300 kg, it is ideal for lifting and transferring a person from or to bed, floor, chair, wheelchair, toilet etc.

Molift Quattro can be delivered in two sizes. The length of the posts are fixed. The lowest allowed distance from floor to ceiling is 2400 mm.


Molift Air or Molift Nomad with trolleys are ordered separately.

Accessory



Adjustable Foot

The Molift Quattro can be optimized with a height adjustable foot on each leg..
Item no 1610200

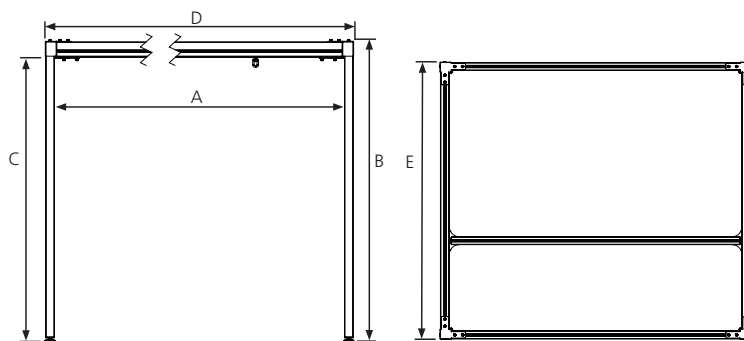

A: Travers lengths	3000 or 4000 mm
B: Height without feet	2330 mm
B: Height - with adjustable feet (accessory)	2380 - 2430 mm
C: Internal height without feet	2190 mm
D: Width	3160 - 4160 mm
E: Depth	3160 - 4160 mm
A x C: Max internal space - without feet	3000 x 2190 mm
A x C: Max internal space - without feet	4000 x 2190 mm
A x C: Max internal space - with adjustable feet	3000 x 2240 m - 2290 mm
A x C: Max internal space - with adjustable feet	4000 x 2240 m - 2290 mm

Product name	Internal space	Item no.	Facts
Molift Quattro 3	3000 x 3000 mm	1600030	Safe Working Load: 300 kg Material: Steel, Anodized Aluminum, Plastic Expected Lifetime: 10 years
Molift Quattro 4	4000 x 4000 mm	1600035	

Weight of unit:

- Skinner Kg/m: 4.2 kg/m/8.8 lbs/m
- Leg with adjustable feet: 8.5 kg/18.7 lbs
- Corner: 2.7kg/6 lbs
- Total weight (Complete unit) (3000 x 3000 mm): 110.95 kg
- Total weight (Complete unit) (4000 x 4000 mm): 131.55 kg

Cleaning

Clean/disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.

For the latest news and continuously updated product information – visit: www.etac.com