



Molift Quick Raiser 205

A sit to stand hoist with
outstanding lifting capacity

 **etac**[®]

Molift Quick Raiser 205

A sit to stand hoist with outstanding lifting capacity

Stable and compact

Molift Quick Raiser 205 is a state of the art sit-to-stand hoist, with an electrically adjustable leg base and an outstanding lifting capacity. The standing hoist is comfortable, stable and safe with exceptional manoeuvrability.

Quick Raiser 205 is compact, stable and exceptionally easy to manoeuvre. The unique wheel position, angled push bar and the extraordinary weight distribution of the base provides a small turning radius and easy manoeuvrability even in very narrow spaces. With the small castor option (50 mm/1.97 inch), it is only 75 mm (2,95 inch) high from the floor, and fits under most beds and chairs while the electrically adjustable legs enables accessibility almost everywhere. The fixed-leg option benefits a very small footprint.

Natural pattern of movement

The Active Lifting Arm (four point sling bar) and the RgoSling Active provides extra support when hoisting from sit to stand. The inclined column encourages a natural pattern of movement and fulfils the hoisting needs of users with some weight bearing capacity. The Active Lifting Arm combined with the sling provides a comfortable, stable and safe hoisting solution.



Molift QuickRaiser 205

Battery	Variant	Adjustable leg	Fixed leg (FL)
SLA battery	Active, small casters	29000	29150
	StandUp, small casters	29001	29151
	Active, large casters	29100	29200
	StandUp	29101	29201
LI battery	Active, small casters	29002	29152
	StandUp, small casters	29003	29153
	Active, large casters	29102	29202
	StandUp, large casters	29103	29203

Model item no.

Technical

Battery SLA 24V	2,9 Ah
Charge time	6 hours
Rated performance	100 hoists/75 kg (100 hoists/165 lbs)
Lifting Range	600 mm (23.62 inch)
Turning radius	1150 mm (45.28 inch)
Hoisting speed 75 kg (165 lbs)	35 mm/s (1.38 inch/s)
Hoisting speed 205 kg (451 lbs)	30 mm/s (1.18 inch/s)

For complete variant matrix, contact your local distributor.



Weight

Maximum user weight (SWL)	205 kg (451 lbs)
Total weight Adjustable leg	44 kg (97 lbs)
Total weight Fixed leg	42,5 kg (93,6 lbs)

Clean the product with a solvent-free cleaning agent with a pH level between 5-9, or with a 70% disinfectant solution.

Dimensions


External width of wheelbase	620 mm (24,4 inch)
Internal width of wheelbase	200 mm (7.87 inch)
Height of footboard:	
Castor size Ø 50/75 mm	65 mm (2,56 inch)
Castor size Ø 75/100 mm	90 mm (3,54 inch)

Features



Ergonomic design

The slanting column lifts the user up and forwards to imitate a natural movement pattern.

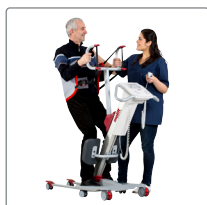

Comfort, safety and stability

The Active Lifting Arm (four point sling bar) and the RgoSling Active provides extra support when hoisting from sit to stand.



Automatic warning

Automatic warning for (annual) periodic inspection and service.



Easy to manoeuvre

Easy to operate in moving and handling situations; toilet, sit-to-stand and stand-up training.



The image includes Active lifting arm & RgoSling Active.

Accessories



StandUp lifting arm

Ideal to use with RgoSling StandUp slings.

StandUp Lifting Arm 2910061



Footplate

This accessory benefits a 10 cm (4inch) elevated footrest. The attachment key is easy accessible and makes it convenient to lock or unlock the accessory to the Quick Raiser.


Footplate 2910065



Set of inlay for footplate

Inlays with slip-resistant material to increase friction.

Set of inlay for footplate..... 2920151

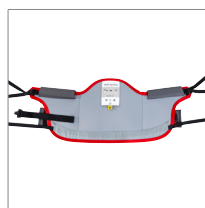

RgoSling Active & RgoSling StandUp

The slings are available in size XS to XL
For item numbers and measurements see www.etac.com



Complete kit Bracket & Cushions

The alternative leg supports benefits an option of small or large cushions, and a leg strap for safety measurements

Kit with bracket & small cushions ... 2910059
Kit with bracket & large cushions 2910056



Set of Cushions

Only contains cushions, no leg straps.

Set of small cushions..... 2910076
Set of large cushions 2910077



Etac is a world-leading developer and provider of ergonomic assistive devices and patient handling equipment. Our heart lies in the solutions that optimise quality of life for the individual, their family and caregivers.

For the latest news and continuously updated product information – visit: www.etac.com

Etac AB
Kista Science Tower
SE-164 51 Kista, Sweden
Tel +46 371 - 58 73 30
info@etac.se www.etac.com

 **etac**[®]
Creating Possibilities